



DRAWING NO. CV15700LAPX	REV G	Customer Drawing	REVISION HISTORY			
			REV	DESCRIPTION (ECO)	DATE	APPROVED
			F	ADD SUFFIX "GU" GREEN URETHANE ECO 6516	5/19/15	AK
			G	REMOVE SUFFIX "C" CLEAR PVC ECO 6576	7/14/15	AK

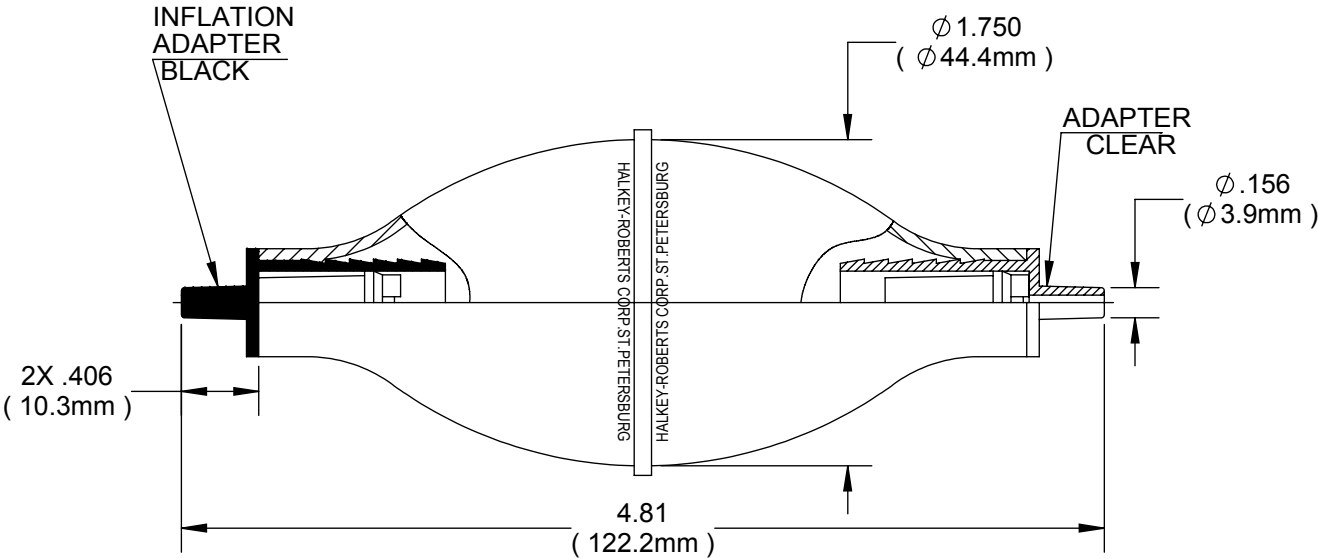


NOTES:

- SUFFIX "X" INDICATES: HAND PUMP COLOR & MATERIAL.
DB - DARK BLUE, PVC
CU - CLEAR, URETHANE
GU - GREEN URETHANE
DBU - DARK BLUE, URETHANE
DBU2 - DARK BLUE, URETHANE WITH 12.5 PSIG TEST AS FOLLOWS.
- ASSEMBLY MUST BE CAPABLE OF DELIVERING 12 PSIG MINIMUM.
TEST PRESSURE: 12.5 PSI
TEST TIME: 5 SECONDS
MAX PRESSURE DROP: .5 PSI.
- CRITICAL DIMENSIONS ARE CIRCLED AND SUBJECT TO INSPECTION.
- DIMENSIONS IN PARENTHESES ARE FOR REFERENCE ONLY.
- ALL DIMENSIONS ARE IN INCHES, WITH NOMINAL IN METRIC BELOW.

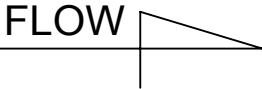

MOLDED PARTS:
GATES, PARTING LINES, AND FLASH PER QUALITY CRITERIA

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS SOLE PROPERTY OF HALKEY-ROBERTS CORPORATION ST. PETERSBURG, FL. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF HALKEY-ROBERTS CORPORATION IS PROHIBITED.

DRAWN BY	M.S.	DATE	04-13-89
APPROVED BY		DATE	
MATERIAL	SEE NOTES		
FINISH			
UNLESS OTHERWISE SPECIFIED:		THIRD ANGLE PROJECTION	
DIMENSIONS ARE IN INCHES			
TOLERANCES:			
ANGULAR: ± 1°			
TWO PLACE DECIMAL ± .01			
THREE PLACE DECIMAL ± .005			

Halkey-Roberts

INFLATE, DEFLATE PUMP

SIZE	DRAWING NO.	REV
A	CV15700LAPX	G

SCALE	DO NOT SCALE DRAWING	SHEET
1:1		1 OF 2


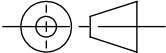
DRAWING NO. CV15700LAPX	REV G	Customer Drawing	REVISION HISTORY			
			REV	DESCRIPTION (ECO)	DATE	APPROVED
				SEE SHEET ONE		

BILL OF MATERIAL:

INFLATION ADAPTER: BLACK PVC , OR BLACK POLYCARBONATE
CATHETER VALVE BODY: POLYPROPYLENE
STEM: SILICONE

HAND PUMP :
"DB" DARK BLUE: 3259X60DB PVC DARK BLUE
"CU" CLEAR: RxT70A URETHANE
"GU" GREEN: RxT70AGU URETHANE
"DBU" DARK BLUE: RxT70ADBM URETHANE
"DBU2" DARK BLUE: RxT70ADBM URETHANE

ADAPTER: CLEAR PVC, OR CLEAR POLYCARBONATE
CATHETER VALVE BODY: POLYPROPYLENE
STEM: SILICONE

MOLDED PARTS: GATES, PARTING LINES, AND FLASH PER QUALITY CRITERIA PROPRIETARY AND CONFIDENTIAL THE INFORMATION CONTAINED IN THIS DRAWING IS SOLE PROPERTY OF HALKEY-ROBERTS CORPORATION ST. PETERSBURG, FL. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF HALKEY-ROBERTS CORPORATION IS PROHIBITED.	DRAWN BY M.S.		DATE 4/13/89	 Halkey-Roberts	
	APPROVED BY		DATE		
	MATERIAL SEE BOM				
	FINISH			TITLE INFLATION, DEFLATION PUMP	
UNLESS OTHERWISE SPECIFIED: DIMENSIONS ARE IN INCHES TOLERANCES: ANGULAR: ± 1° TWO PLACE DECIMAL ± .01 THREE PLACE DECIMAL ± .005	THIRD ANGLE PROJECTION		SIZE A	DRAWING NO. CV15700LAPX	REV G
			SCALE 1:1	DO NOT SCALE DRAWING	
				SHEET 2 OF 2	